


Bitvae

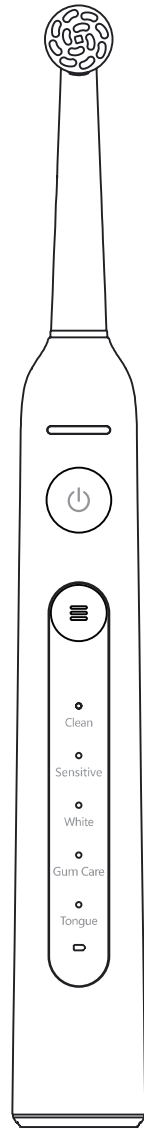

**Daily R2 Rotary
Electric Toothbrush**

SERVICE SUPPORT

1. A 1-year warranty begins on the day of purchase and covers defects caused by non-human factors.
2. Within the warranty period, we will replace the parts or the whole product in case of any defect arising from materials or product malfunction.
3. The following events are not covered by the warranty: brush head defects; damage caused by improper use, for example, improper use of the power supply, plugged into an improper socket, artificial damage, etc.; damage due to the use of the product for other purposes or the use of components not supplied by us; aging and wear caused by normal use. If you have any questions, please contact us at any time.

ENVIRONMENTAL PROTECTION



There is a built-in rechargeable lithium battery in the product. When disposing of the toothbrush, do not mix it with general household waste. It should be sent to the battery recycling center for disposal in accordance with local disposal management regulations.

The symbol on the product or on its packaging indicates that this product may not be treated as household waste. Instead it should be handed over to the applicable collection point for

the recycling of electrical and electronic equipment.

By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. For more detailed information about the recycling of this product, please contact your local council office or your household waste disposal service.

CONTACTS

After-sale Guarantee

Any questions, please contact us. We will provide a solution within 24 hours.

✉ support@bitvae.com

☎ +1 972 679 9068

🌐 <https://bitvae.com>

Clevo Innovation INC

ADD: 2060 N COLLINS BLVD STE 110 RICHARDSON TX 75080

MADE IN CHINA

UK REP

APEX CE SPECIALISTS LIMITED

89 Princess Street, Manchester, M1 4HT, UK

EC REP

APEX ES SPECIALISTS, S.L.

Calle Puerto de la Morcuera, 13, PL. 1 OF. 18-8 28919

Leganés Madrid



Points de collecte sur www.quefairemesdechets.fr
Privilégiez la réparation ou le don de votre appareil !

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING (THIS APPLIANCE)

1. For the first use, it is recommended to start the toothbrush in the 'Sensitive' mode position. If you experience gum discomfort or bleeding gums, please stop using the product and consult your dentist.
2. For the first use, do not put the toothbrush in your mouth before confirming the battery levels. Turn on the toothbrush to make sure the product is fully charged. If the toothbrush has a low battery or no power at all, please fully charge it before use.
3. Do not use, store or charge the toothbrush next to a fire area or high-temperature atmosphere.
4. Do not place the product in a humid environment or in liquid for a long time to avoid product malfunction.
5. Do not disassemble, repair or alter the product or install other accessories or power sources without authorization and guidance.
6. Prevent the product from falling or strong impact.
7. This product is for oral cleaning only, do not use it for other purposes.
8. Do not use this product while it is charging for safety reasons.
9. Make sure to charge the product with a DC5V, 1A adapter.
10. Make sure the charging port or power adapter is dry enough before charging the toothbrush or connecting it to the power source.
11. Periodically check if the power cord is damaged. If the cable is damaged or the toothbrush does not work, please stop using it.
12. In case of skin contact or damage from battery leakage, please stop using the product immediately.
13. Do not put the product in your mouth or brush your teeth while lying down. This can easily lead to accidental swallowing of parts.
14. Children should not play or disassemble this product. Do not place this product where children or babies can reach it.
15. Small parts such as "brush head protective cap" are included in the product package. Please keep it away from children to prevent them from swallowing these parts.
16. People who have severe periodontal disease or who have undergone oral surgery should consult a doctor before use.
17. Users experiencing problems with pacemakers or other implanted devices should consult a medical professional before using this product.
18. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
19. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
20. This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of.
21. The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
22. The power supply with output power not exceeding 15W in normal and abnormal working conditions must be used.

23. Do not expose the appliance or battery to excessive temperatures.
24. Be aware of the risk of terminals of the battery-operated appliance or battery being short circuited by metal objects.

DANGER- To reduce the risk of electrocution

1. Always charge the product in a safe place away from water.
2. Do not reach for the product that has fallen into the water while charging. Unplug immediately.
3. Never charge with a damaged power cord or plug.
4. Make sure the charging port is dry enough before connecting it to the power supply.


WARNING - To reduce the risk of burns, electrocution, fire, or injury to persons.

1. A product should never be left unattended when plugged in.
2. Close supervision is necessary when this product is used by, on, or near children or invalids.
3. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center for examination and repair.
5. Keep the cord away from heated surfaces.
6. Never use while sleeping or drowsy.
7. Never drop or insert any object into any opening or hose.

SAVE THESE INSTRUCTIONS

Tip: The product, accessories, user interface and other illustrations hereby are schematic diagrams for reference only. There may be a slight difference between the product and the schematic diagram due to the update or upgrade of the product. Please take the actual object as the standard.

1. Product Design Descriptions

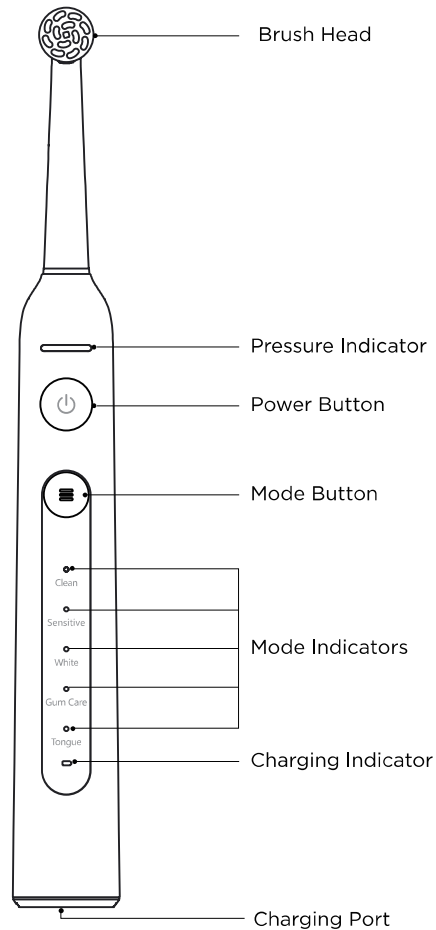

2. Specifications

Product Name	Daily R2 Rotary Electric Toothbrush
Model	R2
Rated Power	2W
Mode	5 Modes (Clean, Sensitive, White, Gum Care, Tongue)
Waterproof Rating	IPX7
Type of Charging	Type-C Direct Current
Input Voltage	5V $\overline{=}$ 1A
Product Measurement	37*29.6*247.5mm
Net Weight	146g
Storage Conditions	Room temperature, dry, ventilated

3. Features

5 modes

Clean: Standard cleaning mode that can effectively remove plaque. Suitable for users in normal oral conditions to maintain healthy oral hygiene.

Sensitive: Sensitive care mode protects the user's gums. Gently and thoroughly cleanses the sensitive areas of the mouth.

White: Whitening and polishing mode that can effectively remove tooth stains and make users' teeth smooth and white.

Gum care: On the premise of optimal cleaning result, the gum care mode massages the gums through the rotation of the frequency conversion toothbrush. This helps stimulate blood circulation which has a positive impact on the gums.

Tongue: Effectively removes bacteria on the tongue coating and improves the overall oral health of the mouth.

Smart pressure indicator

When too much pressure is applied to the teeth, the pressure indicator will activate. It will light up to remind you to reduce pressure on your teeth and prevent damage to your gums.

Vibration gradually increases

To prevent toothpaste splashing, the rotation frequency gradually increases when each mode is activated.

Mode memory function

The default startup mode remains the same as the mode that was set at the time the toothbrush was last powered off.

Smart brushing timer

'Clean', 'Sensitive', 'White' and 'Gum Care' modes have a 3-minute smart timer with 30-second pauses to remind you to change brushing area. The 'Tongue' cleaning mode works for 20 seconds ensuring that you have sufficient time to clean your tongue and remove the white coating from your tongue.

A low-power reminding function

When the battery is running low and you turn on the toothbrush, the charging indicator light will turn red or the toothbrush will not turn on normally as an indicator that the battery is running low and a reminder to charge the toothbrush.

Over-voltage protection

When the input voltage is too high, the toothbrush will enter a self-protection state and the charging process will stop automatically.

Work time

The toothbrush's runtime per charge is typically more than 30 days when fully charged (in "Sensitive" mode for 2 minutes at a time, twice a day).

Charging time

Type C charging, the toothbrush will be fully charged in 3 hours.

Charging indicator

During the charging process, the charging indicator light will flash. When the toothbrush is fully charged, the charging indicator will flash in white and then turn off.

Waterproof

Waterproof design, IPX7 waterproof level.

4. Brushing Guidelines

Charging Instructions

Indicators	Descriptions
The indicator flashes red.	The battery is running low, please charge the toothbrush.
The indicator flashes slowly.	The toothbrush is in the process of charging.
The indicator flashes white and then turns off.	The battery level of the toothbrush is sufficient.

Note: Output cannot exceed 5.5V, otherwise the device will burn out, resulting in a permanent end to its operational life.

Installation or Removal of Brush Heads



Tips: When using the toothbrush, please choose a matching brush head and replace it with a new brush head every 3 months for an optimal oral care experience.

Install the brush head: Align the bottom of the brush head with the shaft, push it in slightly until the brush head is locked in the handle.


Detach the brush head: Hold the toothbrush handle and pull out the brush head along the shaft.

1. Damp the bristles and apply 1 cm of toothpaste to the bristles.


Tips: To prevent toothpaste from falling and splashing, please put the brush head with toothpaste in your mouth before turning on the toothbrush.

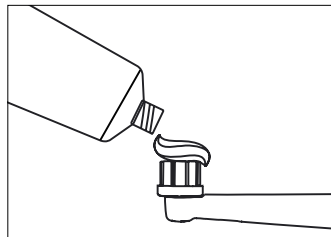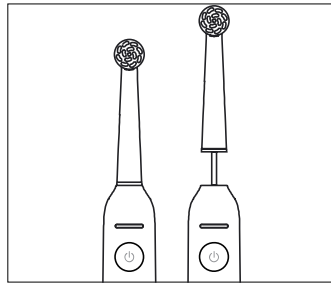


2. Position the bristles towards the outer side of a tooth at a 45° position as shown in the picture, so that the bristles can cover the tooth and the gingival sulcus.

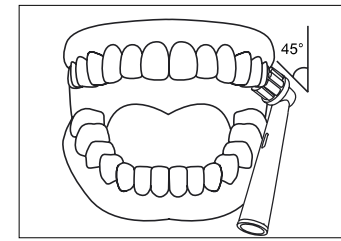

3. Briefly press the power button of the toothbrush, select the appropriate mode, and gently move the bristles from the inner to the outer sides of your teeth.

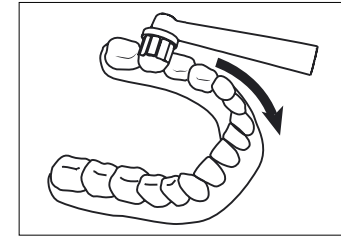

4. Clean each side of the tooth according to the 4 areas shown in the image.

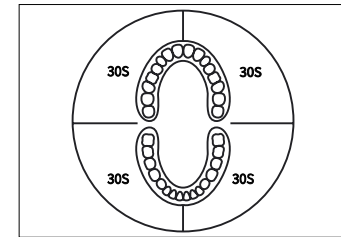

5. After use, rinse your mouth and toothbrush with clean water.

